



CONTROL	L	N	H
---------	---	---	---

IMPORTANT: The barcode is for use only on the CELL-DYN Ruby. Refer to the appropriate System Operator's Manual for proper use of CELL-DYN Calibrator and Control Products.


IMPORTANT: Mixing and Handling

1. Remove a vial of the control from the refrigerator and warm to room temperature (18° to 30° C) for 15 minutes before use.


2. To mix: (**Do NOT mix mechanically or vortex.**)

For a video demonstration, visit www.corelaboratory.abbott and navigate to the Customer Portal → Technical Library → Other Reference Documents → Hematology Aids.

a. Hold the vial vertically and roll each vial between the palms of the hands for 15-20 seconds.

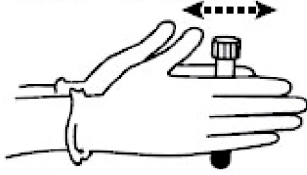

b. Continue to mix by holding the vial by the ends between the thumb and finger, rapidly inverting the vial 20 times end-over-end using a very quick turning motion of the wrist.

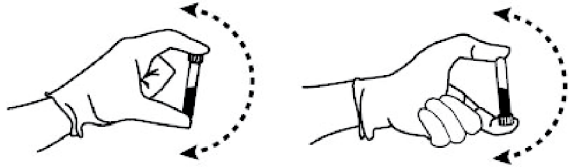

c. Analyze immediately after mixing. Subsequent analyses during this test period may be performed by inverting the vial 5 times prior to instrument analysis.

d. Steps a-c must be repeated upon removing the sample from the refrigerator for the entire open-vial time period regardless of the method of analysis (open tube, cap piercing, auto sample or manual sample).

3. Refer to the appropriate CELL-DYN System Operator's Manual for information about analyzing control specimens.

4. FOR AUTOMATED SAMPLING OR MANUAL CLOSED SAMPLING (CS):

- Refer to the appropriate CELL-DYN Operator's Manual. Remove the vial from the sample handler immediately after sampling.

FOR OPEN-VIAL SAMPLING:

- Aspirate a sample from the vial.
- Carefully wipe the vial rim and cap with a lint-free tissue.
- Replace the cap, ensuring it is on tight.

After sampling, return vial to refrigerator for maximum open-vial stability. If run in the open mode, wipe the threads of both vial and cap before replacing cap and returning to refrigerator.



2025-03-07

8 Consecutive Day Open-Vial Stability

		CONTROL L		CONTROL N		CONTROL H	
		LOT L4358		LOT N4358		LOT H4358	
SYSTEM	PARAMETER	ASSAY VALUE	± MEAN RANGE *	ASSAY VALUE	± MEAN RANGE *	ASSAY VALUE	± MEAN RANGE *
CELL-DYN 3700 SYSTEM	WOC 10 ⁹ /L	N/A	N/A	N/A	N/A	N/A	N/A
	WIC 10 ⁹ /L	N/A	N/A	N/A	N/A	N/A	N/A
	WBC 10 ⁹ /L	N/A	N/A	N/A	N/A	N/A	N/A
	NEU 10 ⁹ /L	N/A	N/A	N/A	N/A	N/A	N/A
	NEU %	N/A	N/A	N/A	N/A	N/A	N/A
	LYM 10 ⁹ /L	N/A	N/A	N/A	N/A	N/A	N/A
	LYM %	N/A	N/A	N/A	N/A	N/A	N/A
	MONO 10 ⁹ /L	N/A	N/A	N/A	N/A	N/A	N/A
	MONO %	N/A	N/A	N/A	N/A	N/A	N/A
	EOS 10 ⁹ /L	N/A	N/A	N/A	N/A	N/A	N/A
	EOS %	N/A	N/A	N/A	N/A	N/A	N/A
	BASO 10 ⁹ /L	N/A	N/A	N/A	N/A	N/A	N/A
	BASO %	N/A	N/A	N/A	N/A	N/A	N/A
	RBC 10 ¹² /L	N/A	N/A	N/A	N/A	N/A	N/A
	HGB g/dL	N/A	N/A	N/A	N/A	N/A	N/A
	HCT %	N/A	N/A	N/A	N/A	N/A	N/A
	MCV fL	N/A	N/A	N/A	N/A	N/A	N/A
MCH pg	N/A	N/A	N/A	N/A	N/A	N/A	
MCHC g/dL	N/A	N/A	N/A	N/A	N/A	N/A	
RDW %	N/A	N/A	N/A	N/A	N/A	N/A	
PLT 10 ⁹ /L	N/A	N/A	N/A	N/A	N/A	N/A	
MPV fL	N/A	N/A	N/A	N/A	N/A	N/A	
CELL-DYN Ruby SYSTEM	WBC (WOC) 10 ⁹ /L	3.8	0.4	6.5	0.7	14.8	2.5
	WBC (NOC) 10 ⁹ /L	3.9	0.4	6.6	1.0	15.0	2.5
	NEU 10 ⁹ /L	2.2	0.3	3.8	0.8	8.6	2.0
	NEU %	58.6	8.0	58.8	8.0	58.4	10.0
	LYM 10 ⁹ /L	1.0	0.3	1.7	0.8	3.8	2.0
	LYM %	25.3	9.0	25.6	9.0	25.9	10.0
	MONO 10 ⁹ /L	0.4	0.2	0.6	0.4	1.5	0.6
	MONO %	9.8	5.0	9.5	5.0	9.9	3.0
	EOS 10 ⁹ /L	0.1	0.1	0.2	0.2	0.4	0.2
	EOS %	3.0	3.0	2.6	2.0	2.6	1.0
	BASO 10 ⁹ /L	0.1	0.1	0.2	0.2	0.6	0.6
	BASO %	3.7	3.0	3.5	3.0	3.3	3.0
	RBC 10 ¹² /L	2.94	0.15	4.21	0.20	5.29	0.30
	HGB g/dL	7.3	0.4	11.3	0.6	15.2	0.8
	HCT %	20.5	1.5	31.5	2.5	41.2	3.5
	MCV fL	69.8	4.0	74.8	4.0	77.8	5.0
	MCH pg	24.8	2.0	26.8	2.0	28.7	2.0
MCHC g/dL	35.6	2.3	35.9	3.0	36.9	2.3	
RDW %	13.8	2.5	12.2	2.5	10.8	2.5	
PLT 10 ⁹ /L	80	20	228	30	527	60	
MPV fL	9.9	3.0	9.7	3.0	9.8	3.0	

* The **MEAN RANGE** does not represent standard deviations (SD).
NOTE: Flags may occur with control materials and should be disregarded.



CELL-DYN and CELL-DYN Ruby are trademarks of Abbott Laboratories in various jurisdictions.

Abbott Laboratories
Diagnostics Division
Abbott Park, IL 60064
USA

EC REP Abbott GmbH & Co. KG
Max-Planck-Ring 2
65205 Wiesbaden
Germany
+49-6122-580

MANUFACTURED FOR
Abbott Laboratories



REF 08H59-01, 08H59-02

